


Aaron® Surgi-Center PRO 200 Watt OR Generator | A2350

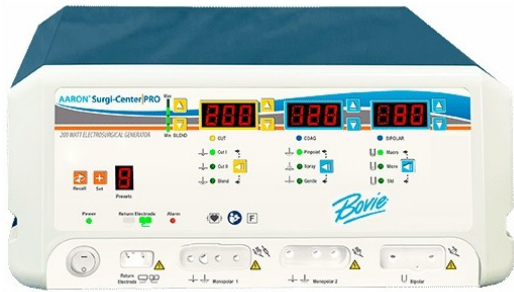

Stock Photo - Options May Vary

FACTORY NEW

The Aaron® Surgi-Center PRO 200 Watt OR Generator is a high-power operating room generator that delivers 12 monopolar and bipolar RF energy modalities for surgical precision. The unit incorporates advanced safety technology including BovieNEM™ return split pad contact quality monitoring and BovieFDFS™ tissue sensing technology that measures tissue 5 million times per second and adjusts to varying impedances 5000 times per second.

The system features variable power adjustment during procedures, simultaneous dual monopolar activation with audible tones, and the ability to store ten user-defined RF presets for easy recall. Bovie's Cut II mode is specifically tailored for laparoscopic procedures with precisely managed cutting energy to minimize the risk of capacitive coupling and alternate site burns.

The unit includes self-test circuits upon power-up and is supported by Bovie's 4-year standard warranty. The compact design weighs less than 20 pounds and measures 37.5 cm wide by 16.5 cm high by 46 cm deep.

FEATURES

- Bovie's Cut II is explicitly tailored for Laparoscopic procedures. The mode offers a precisely managed cutting energy that limits the output, thereby minimizing the potential risk of capacitive coupling and alternate site burns.
- The Aaron SurgiCenter|PRO offers simultaneous activation in spray with audible activation tones. The simultaneous delivery produces equal and even power outputs in each accessory for safer delivery of dual monopolar energy.
- The power can be adjusted up or down in variable steps while performing procedures to increase procedural efficiency and provide the appropriate amount of power to the surgical site.
- The Aaron® SurgiCenter|PRO recalls the last used power and mode. The surgeon can store ten user-defined RF presets for easy recall.
- This return electrode contact quality monitoring system is designed to minimize the risk of patient burns at the return electrode site. This feature works in conjunction with disposable split-return electrodes.
- Bovie FDFS™ - The tissue sensing technology measures tissue 5 million times a second and adjusts to varying impedances 5000 times a second for accurate energy delivery.
- The automatic system self-checks upon powering up the

SPECIFICATIONS

- Audible/Visible Output Meter: Yes
- Pre-set Feature: Yes
- Monopolar HP Ports: 2
- Dual Activation Tone Alert: Yes
- Weight
- Dimensions: W - 37.5 cm (14.75 in) X H - 16.5 cm (6.5 in) X D - 46 cm (18.1 in)
- Warranty: 4 years, Additional Warranty Available

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free 800-627-3215 | Local 813-740-8640

salesteam@venturemedical.com