


Smoke Shark III Smoke Evacuator Unit 120V | SE03

FACTORY NEW

The Aaron Bovie Smoke Shark III is a compact smoke evacuator designed for medical environments requiring efficient smoke removal during electrosurgical procedures. It effectively reduces biological, physical, and chemical hazards while fitting easily into physician offices and mounting conveniently on walls. The unit is engineered for procedures lasting fifteen minutes or less.

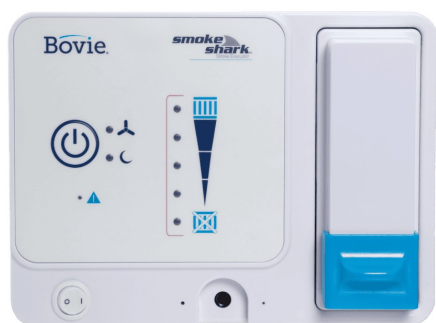

Stock Photo - Options May Vary

The Smoke Shark III features a 4-stage ULPA filtration system that captures 99.999% of particulates and microorganisms as small as 0.1–0.2 microns, maintaining clean air quality and protecting medical staff. The quiet motor enhances patient comfort while maintaining high performance throughout procedures.

The unit includes the SharkSkin Adapter (SS95) for quick setup, the Smart Shark Remote Activation (SERS2) to optimize filter use and extend filter life, and a 6-hour filter (SF06) for high-efficiency filtration in busy medical settings.

FEATURES

- Compact for the physician office
- Quiet and Effective Smoke Evacuation
- Multi-procedure 6-hour filters
- Includes SharkSkin adapter, Smart Shark remote activator, and one filter
- Excellent particulate capture

SPECIFICATIONS

- Height: 7.1" (18 cm)
- Width: 9.4" (24 cm)
- Depth: 6.1" (15.5 cm)
- Weight: 6.25 lb (2.8 kg)
- Line Voltage: 100 to 120 VAC
- Line Frequency: 50 to 60 Hz

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free | Local
800-627-3215 | 813-740-8640

salesteam@venturemedical.com

[View this product online →](#)