


Zoll AED Plus Biphasic AED – New

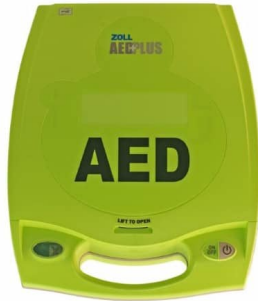

Stock Photo – Options May Vary

FACTORY NEW

The Zoll AED Plus Biphasic is an automated external defibrillator designed for use by both novice rescuers and medical professionals. It provides real-time CPR coaching and feedback, monitoring compression rate and depth to guide rescuers through proper technique. Audio and visual prompts deliver immediate guidance, including alerts to push harder when compressions are insufficient and confirmation of "Good compressions" when performed correctly.

The device features a single-piece electrode pad for rapid and accurate placement, with a 5-year shelf life. It operates on consumer-grade lithium batteries available from retail stores and is highly water-resistant, making it suitable for various environments.

The AED Plus uses Zoll's Rectilinear Biphasic waveform technology and includes Intelligent Pediatric Capability to support treatment of pediatric patients. It automatically selects appropriate energy levels and charges in less than 10 seconds, enabling fast response to sudden cardiac arrest.

FEATURES

- Real CPR Help for rate and depth of compressions
- Support for the complete Chain of Survival
- Help to all victims of SCA, even those for which no shock is advised
- A one-piece electrode pad for fast and accurate placement
- Uses lithium camera batteries available from retail stores
- Intelligent Pediatric Capability
- Rectilinear Biphasic Waveform
- Able to follow Guidelines 2010

SPECIFICATIONS

- Dimensions: 5.25"h x 9.5"w x 11.5"d
- Energy Selection: Automatic preprogrammed selection (120J, 150J, 200J)
- CPR: Metronome Rate: Variable 60 to 100 cpm
- Charge Time: Less than 10 seconds with new batteries
- Charge Hold Time: 30 seconds
- Waveform: Rectilinear Biphasic
- Weight: 6.7 lbs
- Power: 10 Replaceable types 123-A lithium batteries

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free | Local
800-627-3215 | 813-740-8640

salesteam@venturemedical.com

[View this product online →](#)