


Mindray Passport 12 Anesthesia Monitor with AG Module - Refurbished

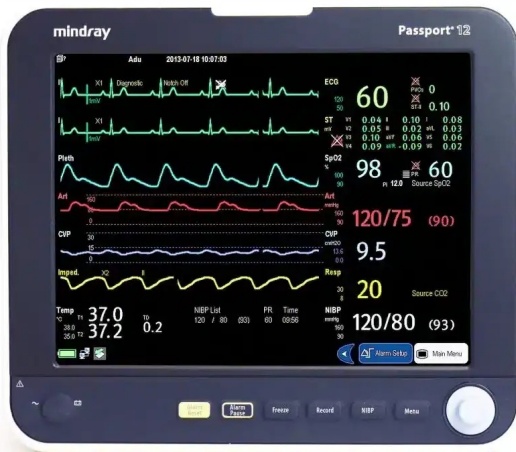

Stock Photo - Options May Vary

FULLY REFURBISHED · CLINICALLY CERTIFIED

The Mindray Passport 12 is a versatile, modular patient monitor designed for acute-care environments including operating rooms and critical care settings. Its plug-and-play architecture supports advanced monitoring capabilities including CO₂, dual invasive pressures, and cardiac output, with an integrated anesthesia gas module for monitoring five anesthetic agents and end-tidal CO₂. The monitor features a 12-inch high-resolution color touchscreen display with up to 8 waveforms and standard parameters including 3- or 5-lead ECG, Masimo SET SpO₂, dual IBP, NIBP, respiration, temperature, and integrated recorder.

The Passport 12 offers extensive data storage for trend data, alarms, events, and NIBP measurements up to 48 hours of full disclosure. It supports both wired and wireless networking options with integrated HL7 and DIAP communication protocols, enabling connectivity to the BeneVision Central Station and Panorama Patient Monitoring System for centralized patient data management. The monitor is compatible with Mindray eGateway for direct HL7 communication to electronic medical records.

Designed for mobility and flexibility, the Passport 12 includes an ergonomic carrying handle and multiple mounting options. Optional features include Nellcor OxiMax SpO₂, Sidestream or Microstream CO₂ monitoring, and ST and arrhythmia analysis capabilities.

FEATURES

- 12" high resolution touchscreen display
- Up to 8 waveforms
- Standard features include 3 or 5-lead ECG, Masimo SET® SpO₂, 2 IBP, NIBP, respiration, temperature, integrated recorder
- Extensive data storage capability for trend data, alarms, events, NIBP measurements and up to 48 hours of full disclosure
- Available wireless 2.4 Ghz (802.11b/g/n)
- Connects to DPM Central Station or Panorama® Patient Monitoring Network
- Integrated HL7 and DIAP protocol interfaces
- Ergonomic carrying handle and multiple mounting options offer unprecedented mobility

SPECIFICATIONS

- Monitor Size: 10.9 H x 5.4 D x 12.7 W in
- Weight: 10 lb (standard parameters configuration, including 2 lithium batteries and a recorder)
- AC Power: 100 to 240 VAC (±10%)
- Current: 1.3 to 0.5 A
- AC Power Input Frequency: 50/60 Hz (±3Hz)
- Temperature Range (Operating): 10°C to 40°C
- Humidity Range (Operating): 15% to 95% non-condensing
- Altitude Range (Operating): (-1,000 to 9,841 feet ASL) | 1,074 hPa to 700 hPa | (806 mmHg to 524 mmHg)

Toll-Free 800-627-3215 | Local 813-740-8640

salesteam@venturemedical.com

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free 800-627-3215 | Local 813-740-8640

salesteam@venturemedical.com