


Pedigo Wall Mounted ivNow®-6 Modular Thermal Fluid Warmer | P-2306-A

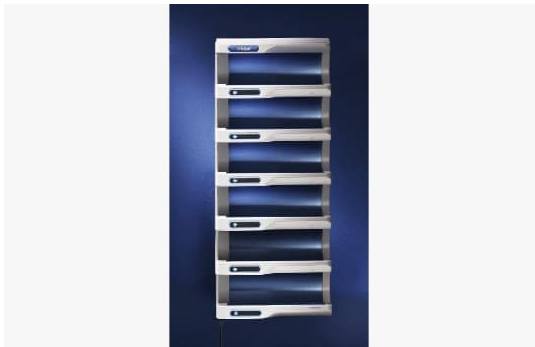

Stock Photo - Options May Vary

FACTORY NEW

The Pedigo Wall Mounted ivNow-6 Thermal Fluid Warmer is a compact, wall-mounted unit designed to warm and maintain safe temperatures of intravenous and irrigation fluids. The warmer features six 1-liter bag capacity with individual control and LED display for each cavity. A sensor in the heating plate detects the presence of a bag and automatically activates the heating mechanism.

Two temperature sensors continuously monitor each fluid bag's temperature and automatically engage the heater as needed to maintain fluid temperature within $+0/-2^{\circ}\text{C}$ ($+0/-3^{\circ}\text{F}$) of the setpoint. This ensures fluids remain at the desired temperature for optimal clinical use without the expense and setup time required by in-line warming systems.

The unit keeps fluids at proper temperature and ensures immediate availability at or near body temperature instead of room temperature. Wall mounting using bracket kits provides convenient placement where needed. This unit is not approved for warming blood or blood products.

FEATURES

- Wall-mounted using Wall Mounting Bracket Kits
- Quickly warms and maintains safe temperatures of intravenous and irrigation fluids
- Specially contoured warming surface cradles one (1) liter solution bag
- An individual control for controlling each cavity, and L.E.D. display
- A sensor in the heating plate detects the presence of a bag and engages the heating mechanism to warm the fluid quickly
- Two temperature sensors continuously read the temperature of the fluid bag
- Sensors auto-engage the heater as necessary to maintain the fluid temperature within $+0/-2^{\circ}\text{C}$ ($+0/-3^{\circ}\text{F}$) of the setpoint temperature
- Pedigo's ivNow fluid warmers are small, convenient units

SPECIFICATIONS

- Dimensions: 35"H x 14.5"W Capacity: 6 bags (liter)
- Display: LED
- Temp Sensors: 2
- Accuracy: Within $+0/-2^{\circ}\text{C}$ ($+0/-3^{\circ}\text{F}$) of the setpoint temperature
- Note: Not approved for blood or blood products

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free | Local
800-627-3215 | 813-740-8640

salesteam@venturemedical.com