


SchureMed Adjustable Foot Board/Table Extension | 800-0285

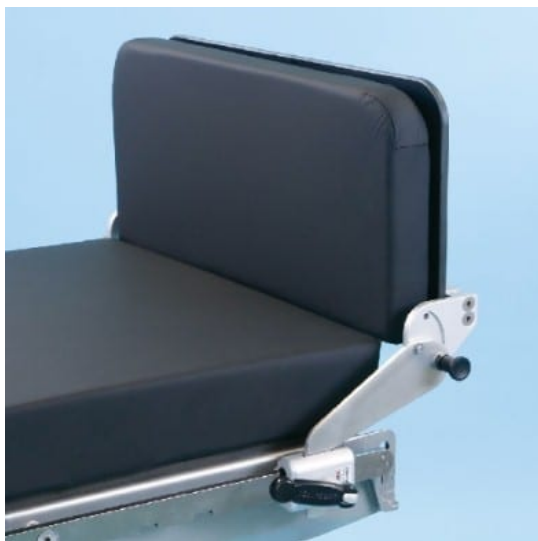

Stock Photo - Options May Vary

FACTORY NEW

The SchureMed Adjustable Foot Board/Table Extension (Model 800-0285) is a versatile surgical table accessory that mounts horizontally as a table extension or vertically as foot support during reverse Trendelenburg positioning. The device features a 5/8" stainless steel mounting post with deluxe foam pads included for patient comfort and safety.

Operation is straightforward: pull both locking knobs to lower the table extension pad to align with the O.R. table pad level, or pull the knobs to raise the extension to a 90° position. The adjustable design accommodates various surgical positioning needs while maintaining stability and ease of use.

The extension supports a 500 lb. proportional patient load and is constructed without natural rubber latex. Installation is designed for single-person setup at the surgical table foot section. The device carries a two-year warranty against manufacturing defects and requires only routine cleaning for maintenance.

FEATURES

- Mounts horizontally as a table extension
- Mounts vertically as a foot support
- Locking knobs to lower table extension pad to the level of the O.R table pad
- Can be raised to a 90° position
- 5/8" Stainless Steel Mounting Post
- Deluxe foam pads are included
- The device supports a 500 lb. proportional patient load (6'4" tall patient per 99% human body model)
- Product is not made with (NRL) Natural Rubber Latex

SPECIFICATIONS

- Length: 24" +/- 0.5"
- Width: 18" +/- 0.5"(Extended)
- Depth: 10" +/- 0.5" (W/Pad)
- Device Weight: 13 +/- 0.5 lbs
- Attaches to the surgical table at the foot section
- Single-person installation

**Brochure and images are a general overview of the product. Actual product options may vary.*

Toll-Free 800-627-3215 | Local 813-740-8640

salesteam@venturemedical.com

[View this product online →](#)