

ECHOCARDIOGRAPHY

Ultrasound Table

- Sliding cardio panel allows you to adjust the opening to the size you need
- Electric or manual sliding cardiac panel opens as wide as 15" (38cm). Great for left or right handed scanning
- 22" to 38" (56-97cm) height range supports ergonomic working heights for both echocardiograms & stress echo
- Designed specifically for vascular & cardiac exams
- Dual tower design for added strength & rigidity allowing a 550 lb (250kg) lift/load capacity



shown in Clay fabric

Sliding Cardio Panel
opens as wide as 15" (38cm)



Inspired by the leading ergonomic experts in sonography

SOUND | ERGONOMICS

- Decrease sonographer injury
- Features minimize exam times



FM584901

made in the USA
with US & imported parts



OAKWORKS®
MEDICAL
Position for Success

ECHOCARDIOGRAPHY

Ultrasound Table

Recommended for:

- Cardiac exams
- Vascular exams
- All abdominal & superficial structures exams
- Musculoskeletal exams
- Biopsies



Adjustable Access/Back Support



ECHOCARDIOGRAPHY TABLE #65258-T*

FEATURES	
2 Section Top	30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range 22" - 38" (56-97cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 15^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Manual Sliding Cardio Panel opens 15" (38cm)
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Adjustable access back support (right) 0-85°, Manual sliding cardio panel (left) opens 15" (38cm), Paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.		#
Accessory Mounting Plate	Fowler end	62123
Central Base Locking System	Not available with Patient Lift Compatible Base	61777
Foot Control	Height, Trendelenburg, Backrest	61901
Foot Control	Height, Trendelenburg, Backrest, Cardio Panel	65679
Non-USA electric system	please contact us for order details	73134
Patient Lift Compatible Base	Not available with Central Base Locking System	65893
Side Rails (set of 2)	23.5" (60cm) long x 6" (15cm) tall (from surface of table)	61781
Sliding Cardio Panel - Electric	Upgrade your manual panel opens 15" (38cm)	73131
ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	62192
T-rail Adapter 8" (20cm) long	Attaches to the backrest for use with USA format T-rail attachable accessories.	61902
Table Extender	24" (61cm) wide x 12" (30cm) long (requires mounting plate)	63487
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support (requires mounting plate)	PKG7554-T*

*Fabric color code must be specified.

ULTRASOUND FEATURES & OPTIONS

CENTRAL BASE

LOCKING SYSTEM #61777

Easily lock your table into place with our unique system that allows locking from both sides of the table.



OPTIONAL

SIDE RAILS #61781

Increase stability and security for your patient. Fold down out of the way when not in use. 23.5" (60cm) long x 6" (15cm) tall (from surface of table)



OPTIONAL

FULL FUNCTION FOOT CONTROLS

#enquire for part number

You can add a foot control that controls all functions.



Multi-function foot control

OPTIONAL

PATIENT LIFT COMPATIBLE BASE #65893

This option allows the use of patient lifts and easy wheel chair transfers to minimize or eliminate caregiver injuries associated with patient transfers. The high clearance base permits all patient lifts to slide easily under the center of the table for stable transfers while not increasing the low height range of the tables that enable easy wheel chair transfers. Must be ordered with table.

- Includes larger, 5" (13cm) diameter Heavy Duty locking casters
- Provides 6.5" (17cm) clearance for patient lift
- Central Base Locking System not available with this option



OPTIONAL

ADJUSTABLE ACCESS BACK SUPPORT

Right: #73132

Adjustable with 5 positions 0-85°, allowing easy access for sonographer to the far side and back support for patients. 10" (25cm) wide x 15.5" (39cm) long



STANDARD

ELECTRIC SLIDING CARDIO PANEL #73131

The electric version allows you to adjust the opening to the size you need, up to 15" (38cm). 10" (25cm) wide x 28" (71cm) long. Great for left or right handed scanning. Can also be used as an armrest.



OPTIONAL

ACCESSORY MOUNTING PLATE

#62123

Add a face rest or table extender to your table.



OPTIONAL

ULTRASOUND ACCESSORIES

DETACHABLE ARM BOARD #62192-T*
Attaches to the T-rail Adaptor and rotates over 140 degrees. 5"(13cm) wide x 20.5" (52cm) long (Requires T-Rail Adaptor)



T-RAIL ADAPTOR #61902
Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) long

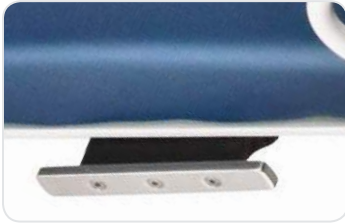


TABLE EXTENDER #63487-T*
Lengthens the table by 12" (31cm) to 84" (213cm) allowing for access to the patient's neck for carotid artery exams and thyroid exams. 24" (61cm) W x 12" (31cm) L



QUICKLOCK™ HEAD REST & AERO-CEL PAD #PKG7554-T*
The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments. (Requires Accessory Mounting Plate)



Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.

PADDING



3" Comfort Foam™
More comfort and support

Pure White-T26	Butter-T09	Saffron-T21	Opal-T12	Clay-T18	Taupe-T28	Earth-T07
Espresso-T27	Ruby-T08	Orchid-T19	Heron-T05	Ocean-T20	Sapphire-T17	Sky Blue-T04
Seafoam-T30	Sage-T22	Blue Grass-T24	Forest-T14	Stone-T13	Pewter-T29	Coal-T01

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

CONTACT:

View our full line of Medical products:

<https://www.oakworksmed.com/pdfdocs/oakworks-medical-catalog/>

717.235.6807 www.oakworksmed.com email: customerservice@oakworks.com