


Ergonomic Advantage Line



General 3 Section EA Ultrasound Table

Our EA Line of tables incorporates all of the ergonomic features and comfort our tables have had for years with electric features that speed up exam times and make positioning patients easier.

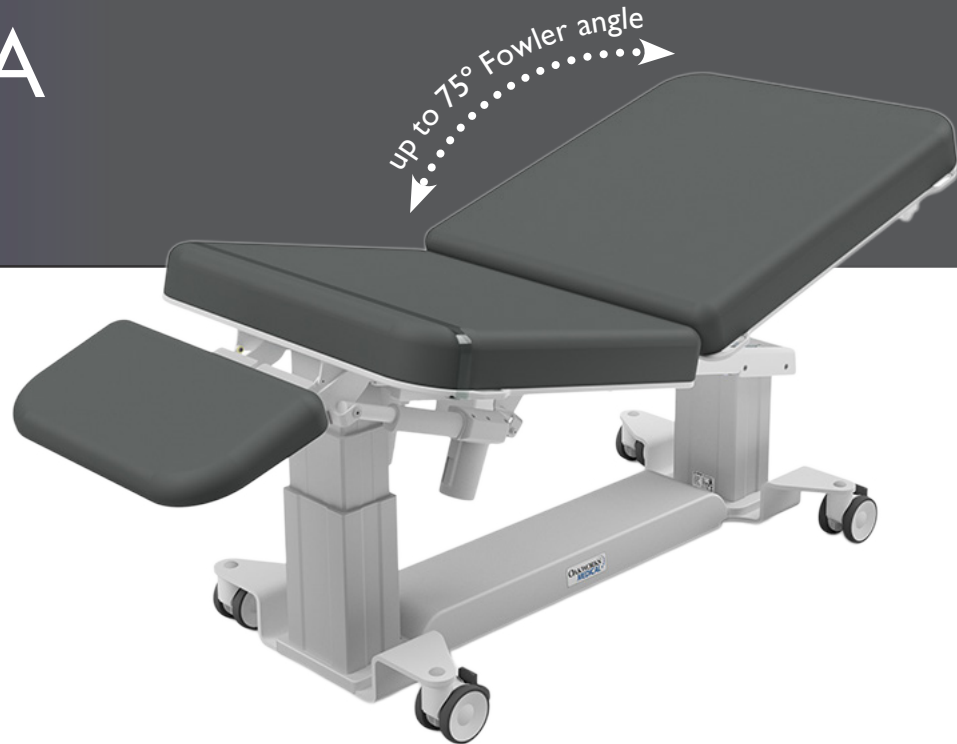


The Oakworks General 3-Section Ultrasound Table offers versatile positioning with electric height, backrest, and leg section adjustments to support a variety of procedures. Its durable design, smooth operation, and patient-friendly features ensure comfort and efficiency in everyday clinical use.

made in the USA
with US & imported parts



Fully Electric Table

General 3 Section EA Ultrasound Table

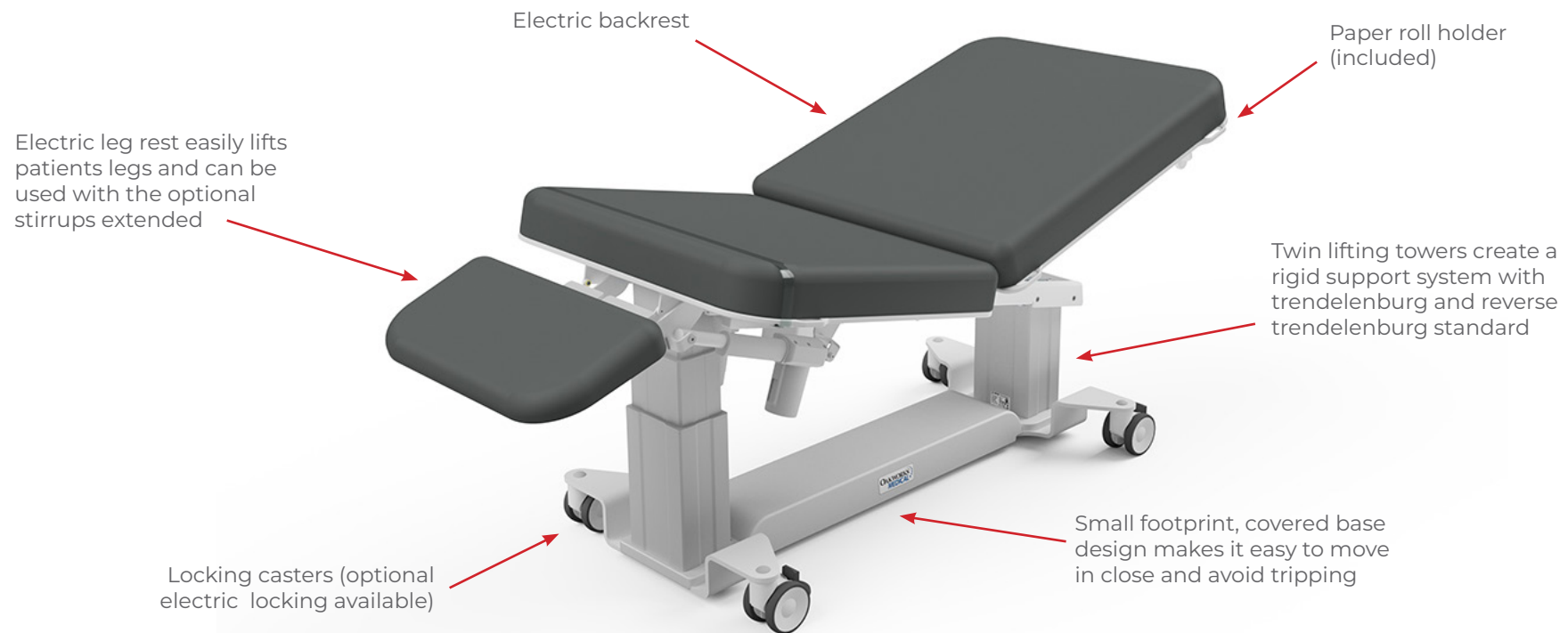

Features

- Designed to use as a table or an exam chair with many options to fit your needs.
- Sleek base design promotes unobstructed compatibility with ultrasound equipment ensuring better ergonomic conditions for the sonographer.
- Designed for versatility—functions as both a table and an exam chair with multiple positioning options.

Recommended for:

- Abdominal & superficial structures exams
- OB/GYN exams (limited use w/ optional stirrups required)
- Vascular exams
- Musculoskeletal exams

Inspired by leading ergonomic experts in sonography

- Decrease sonographer injury
- Minimize exam times

FEATURES	
3 Section Top	26" (66cm) or 30" (76cm) width x 75" (191cm) length
Backrest	34" (86cm) long
Seat	26" (66cm) long
Leg Rest	22" (56cm) width x 11" (28cm) length
Motion 1	Electric Height Range: 19" - 33" (48cm-84cm) low height or 22" - 36" (56-91cm) high height
Motion 2	Electric Trendelenburg / Reverse Tilt $\pm 18^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Electric Leg Rest 0-90°
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Standard TerraTouch™ or Optional Ultraleather
Included Features	4" (10cm) locking casters, paper roll holder
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.

Left Cardio & Right Access / Back Support Panels (available only on 30" tables)	Electric Sliding Cardio Panel opens 15" (38cm) Left side electric access panel opens 10" (23cm) width x 22" (56cm) length
Foot Control	Height, Trendelenburg, Fowler
Powered Central Base Locking	Lock and unlock using hand control
Procedure Tray (available only on 30" tables)	12.25" (6cm) width x 10.25" (26cm) length and holds 2 quarts (1.9 liters)
Retractable footrests	Stirrups
Side Rails (set of 2)	25.5" (65cm) length x 6" (15cm) height (from surface of table)
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.

Detachable Arm Board	5" (13cm) width x 20.5" (52cm) length (Requires T-Rail Adapter)
QuickLock™ Head Rest &	Comfortable, adjustable prone or supine support
Table Extender	24" (61cm) width x 12" (30cm) length
T-Rail Adapter 8" (20cm) long	Attaches to the backrest of the table for use with IV poles, arm rests, and other attachable devices.

*Fabric color code must be specified.

General 3 Section EA Features

ELECTRIC LEG REST



Powered leg rest easily lifts patients legs and can be used with the optional stirrups extended.

BATTERY BACK UP



The battery backup system is capable for running for several hours if the power fails.

HAND CONTROL



Graphic interface hand control with presets quickly positions the patient by controlling multiple functions at once.

REPLACEABLE PADS



Our design includes easily replaceable pads, extending the functional life of the tables while keeping maintenance costs low.

TENTE CASTERS



Rust-proof swivel casters with a total lock mechanism, positioned at the trailing end for enhanced mobility and a high level of cleanability.

ELECTRIC BACKREST



The electronics with graphic interface hand control allow you to easily and quickly change the position of the backrest.

General 3 Section EA Accessories & Options

T-RAIL ADAPTER (ACCESSORY)



Easily installed and removable from multiple dedicated mounting locations on the backrest or seat. Supports up to 150 lbs and fits any 8" (20 cm) American standard accessory.

TABLE EXTENDER (ACCESSORY)



Easily attaches to the head of the table to increase the overall length and accommodate taller patients or specialized positioning needs. Measuring 24" (61 cm) wide by 12" (31 cm) long.

QUICKLOCK™ HEAD REST & AERO·CEL PAD (ACCESSORY)



The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments.

DETACHABLE ARM BOARD (ACCESSORY)



(Requires T-Rail Adapter)

Attaches to the T-rail Adapter and rotates 140 degrees. 5" (13 cm) width x 20.5 (52 cm) length.

ELECTRIC CENTRAL BASE LOCKING (OPTION)



The system allows you to easily lock or unlock the base by simply pushing a button on the hand control.



SIDE RAILS (OPTION)



Our design recesses the side rails under the top when not in use, enhancing ergonomics.

BATTERY BACK UP (OPTION)



A battery back up system capable of running for several hours if the power fails.

ELECTRIC CARDIO & ACCESS PANELS (OPTION)



Electric cardio panel slides open to the size needed for support. Electric access panel offers better access minimizing sonographer reach.

STIRRUPS & PROCEDURE TRAY (OPTIONS)



Adjustable stirrups store neatly under the table. Pair with our removable, retractable tray for a flexible setup suited to various clinical procedures.

AVAILABLE FABRIC COLORS

To view our current list of available fabric colors, scan the QR code or visit:

www.oakworksmmed.com



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible.