

Ergonomic Advantage Line



Vascular with Fowler EA Ultrasound Table

Our EA Line of tables incorporates all of the ergonomic features and comfort our tables have had for years with electric features that speed up exam times and make positioning patients easier.




The Oakworks Vascular Table is built for precision and access, featuring electric height, backrest, and Trendelenburg adjustments to support vascular imaging and procedures. Its ergonomic design and stability enhance workflow efficiency and patient safety.

made in the USA
with US & imported parts



Fully Electric Table

Vascular with Fowler EA Ultrasound Table

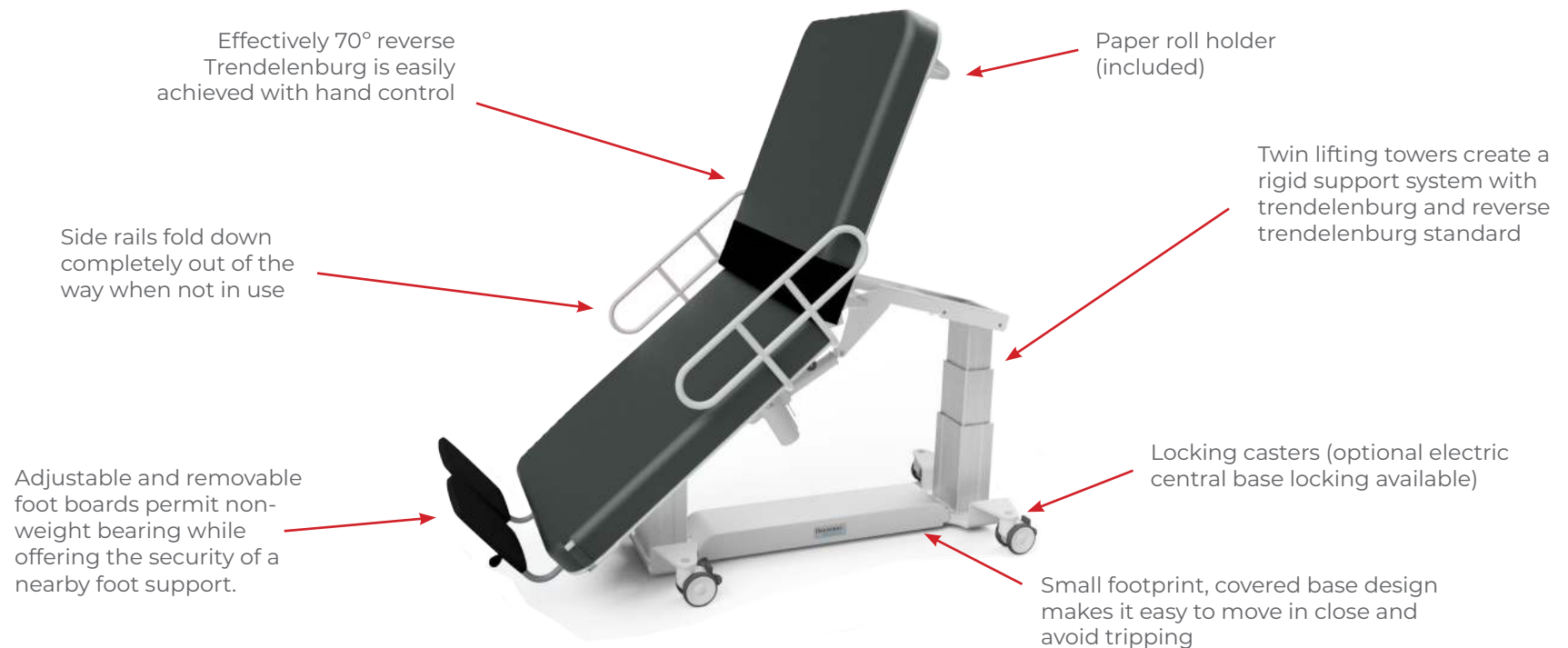

Features

- Safely position clients for venous reflux studies
- Dual towers & electric fowler makes effectively 70° reverse Trendelenburg a reality
- Split foot board design allows non-weight bearing exams and fold down out of the way for carotid artery exams

Recommended for:

- Venous reflux studies
- All other vascular exams
- Abdominal & superficial structures exams
- Musculoskeletal exams
- Biopsies

Inspired by leading ergonomic experts in sonography

- Decrease sonographer injury
- Minimize exam times

FEATURES	
2 Section Top	26" (66cm) or 30" (76cm) width x 72" (183cm) length
Backrest	34" (86cm) long
Seat	38" (97cm) long
Motion 1	Electric Height Range: 19" - 33" (48cm-84cm) low height or 22" - 36" (56-91cm) high height
Motion 2	Electric Trendelenburg 18° Reverse Trendelenburg 40°
Motion 3	Electric Additional Reverse Trendelenburg 0-25° (40° total)
	Electric 0-75° Backrest
Hand Control	All powered motions
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Standard TerraTouch™ or Optional Ultraleather
Included Features	4" (10cm) locking casters, integrated foot board (split, drop away), paper roll holder, patient strap, side rails (2) 25" (64cm) length x 6" (15cm) height (from surface of table)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.

Cardio Panel and Powered Access / Back Support	Electric Sliding Cardio Panel opens 15" (38cm) - left side Electric Access Panel opens 10" width (25cm) x 22" length (56cm)-right side
Powered Central Base Locking	Lock and unlock using hand control
Foot Control	Height, Trendelenburg, Backrest
Battery Backup	Capable of running during power failure
Non-USA Electric System	Please contact us for order details

ACCESSORIES order any time.

Detachable Arm Board	5" (13cm) width x 20.5" (52cm) length (Requires T-Rail Adapter)
Trail Adapter 8" (20cm) long	Attaches to the backrest of the table for use with IV poles, arm rests, and other attachable devices.
Table Extender	24" (61cm) width x 12" (30cm) length
QuickLock™ Head Rest & AeroCel Pad	Comfortable, adjustable prone or supine support

*Fabric color code must be specified.

Vascular with Fowler EA Features

FOOT REST DESIGN



Oakworks exclusive, removable split foot rest design allows you to stagger them, allowing weight or non-weight bearing studies while offering the security of a nearby foot support for the unsupported foot. Foot rests are completely removable and can be folded down out of the way.

HAND CONTROL



Graphic interface hand control with presets quickly positions the patient by controlling multiple functions at once.

REPLACEABLE PADS



Our design includes easily replaceable pads, extending the functional life of the tables while keeping maintenance costs low.

TENTE CASTERS



Rust-proof swivel casters with a total lock mechanism, positioned at the trailing end for enhanced mobility and a high level of cleanability.

SIDE RAILS



Our design recesses the side rails under the top when not in use, enhancing ergonomics.

Vascular with Fowler EA Accessories & Options

T-RAIL ADAPTER (ACCESSORY)



Easily installed and removable from multiple dedicated mounting locations on the backrest or seat. Supports up to 150 lbs and fits any 8" (20 cm) American standard accessory.

TABLE EXTENDER (ACCESSORY)



Easily attaches to either the head or foot of the table to increase the overall length and accommodate taller patients or specialized positioning needs. Measuring 24" (61 cm) wide by 12" (31 cm) long.

QUICKLOCK™ HEAD REST & AERO•CEL PAD (ACCESSORY)



The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments.

DETACHABLE ARM BOARD (ACCESSORY)



(Requires T-Rail Adapter)

Attaches to the T-rail Adapter and rotates 140 degrees. 5" (13 cm) width x 20.5 (52 cm) length.

ELECTRIC CENTRAL BASE LOCKING (OPTION)



The system allows you to easily lock or unlock the base by simply pushing a button on the hand control.

BATTERY BACK UP (OPTION)



A battery back up system capable of running for several hours if the power fails.

ELECTRIC CARDIO PANEL (OPTION)



Electric cardio panel slides open to the size needed for support. Electric access panel offers better access minimizing sonographer reach.

AVAILABLE FABRIC COLORS

To view our current list of available fabric colors, scan the QR code or visit:

www.oakworksmed.com



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible.